A E N O

Quick Start Guide


Toothbrush DB1S / DB2S

(Models: ADB0001S,
ADB0002S)


22.08.2023 | v 2.0.0

ENG.....	7
ARA.....	11
BOS.....	14
BUL.....	18
CES.....	23
DEU.....	27
ELL.....	32
EST.....	36
FRA.....	41
HRV.....	46
HUN.....	50
HYE.....	55
ITA.....	60
KAT.....	65
KAZ.....	70
LAV.....	75
LIT.....	79
NLD.....	84
POL.....	89
POR.....	94
RON.....	99
RUS.....	104
SLK.....	109
SLV.....	114
SPA.....	118
SRP.....	123
UKR.....	128
UZB.....	133
ZHO.....	138


A

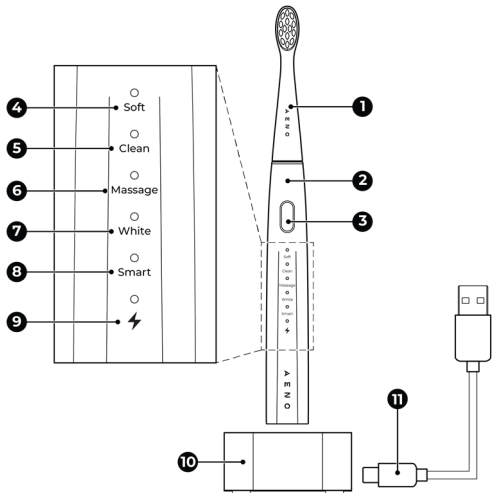
B

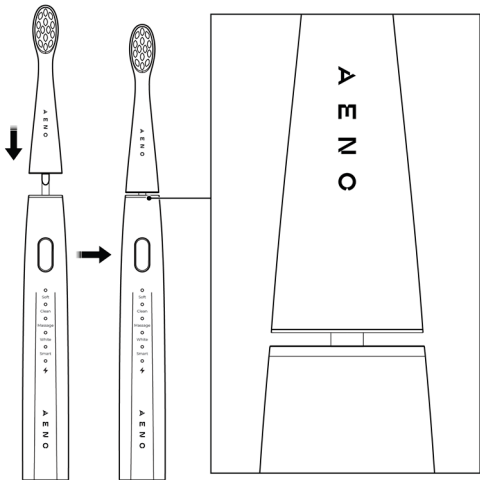C

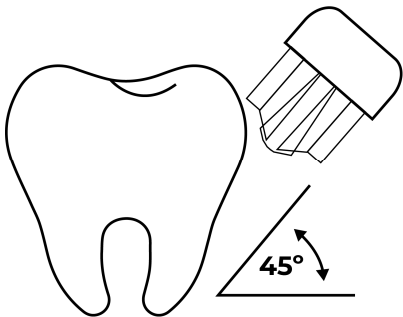D

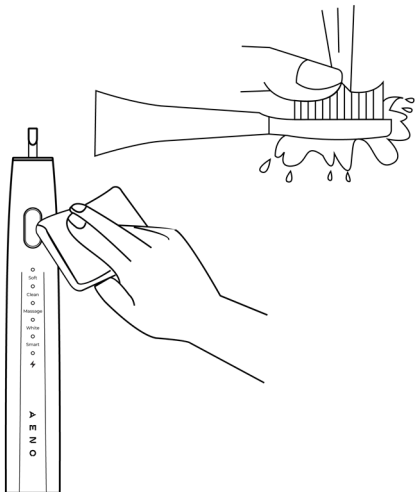ENG The **AENO™ ADB0001S** (white) / **ADB0002S** (black) Smart Toothbrush is designed for daily teeth brushing and gum massage.

Technical specifications

Toothbrush. Battery: non-removable, Li-ion, 1,400 mAh, 3.7 V. Input: 5.0 V (DC), 0.3 A, 1.5 W. Full battery charge time: up to 5.5 h. Operating time after a full charge: up to 90 days. Ingress Protection rating: IPX7. Noise level (at a distance of 0.3 m): ≤ 60 dB. Bristle stiffness: medium. Number of operating modes: 12. Control: manual, remote. BT: 4.2 (2,402–2,480 MHz). Transmit power: +4 dBm / 2.5 mW. Toothbrush size with brush head (diameter × height): 25.8×240 mm. Operating conditions: temperature 0...+40 °C, relative humidity < 90 % (without condensation). Storage conditions: temperature 0...+40 °C, relative humidity < 90 % (without condensation).

Charging base. Input: 5.0 V (DC), 1.0 A, 5.0 W. Ingress Protection rating: IPX5. Size (diameter × height): 45×21 mm.

Scope of supply (see fig. **A**)

Toothbrush (handle (2), 3 replacement brush heads (1), charging base (10)), USB (A) / USB (C) power cable (0.8 m) (11), Quick Start Guide.

WARNING! The power adapter is not included in the package. You can use a power adapter from another device if it provides the required power supply parameters: 5.0 V (DC), 1.0 A.

Safety operation rules

The device is not intended for use by children or by persons with reduced physical, mental or intellectual capabilities if they do not have sufficient experience or knowledge in using the device. If you have a pacemaker, suffer from oral diseases, or if you had oral surgery, consult your doctor before using the device. If your gums begin to bleed after using the device and the bleeding continues for several days, contact your dentist. Do not turn on or use a damaged device. Make sure that the voltage and frequency specified in the technical documentation correspond to the power supply parameters. Charge

If you have any questions or issues with your AENO device, please contact our online support team support@aeno.com or aeno.com/service-and-warranty. This is a better option than visiting the store and will save you time and effort.

the toothbrush only on the original charging base. Avoid dropping the charging base into water. If the charging base accidentally fell into water, disconnect it from the power supply and then remove it from the water and dry thoroughly. Each family member should use an individual replacement brush head. Replace the brush head at least once every three months. Use only the original replacement brush heads. Do not allow children to play with the device and its accessories. When not in use for a long time, charge the device every 2–3 months to prevent a deep discharge of the battery.

Preparing the device for operation

Charging. Unpack the device and place the charging base (10) on a dry, stable horizontal surface. Connect the charging base to a power source using the USB cable (11). Place the handle (2) on the charging base (10). The LED indicator (9) will start to flicker * green.

When the LED indicator (9) turns solid green, the battery is fully charged, the device is ready for operation. When charging is complete, disconnect the charging base from the power supply.

Installing the brush head. Remove the handle (2) from the charging base (10). Place the brush head (1) on the handle (2), pressing lightly until you hear a distinctive click. When properly installed, there should be a small gap between the handle of the toothbrush and the brush head (see fig. B).

LED charge indicator statuses

Indicator colour (9)	Battery charge level
red	low
green	medium/high

Configuring the device operation modes

Turn on the device by pressing the power button (3) once. Subsequent presses of the button will switch operating modes. You can change the mode within 5 seconds after turning on the device. The last used mode is saved in the device settings the next time the device is turned on.

* The indicator is on for 3,000 ms and then is off for 300 ms.

Operating instructions

Wet the bristles of the brush head and apply the desired amount of a dental care product to the bristles. Point the bristles of the brush head towards the surface of the teeth at approximately a 45° angle (see fig. C). Press the power button (3) to turn on the device. Brush your teeth without putting excessive pressure on the enamel. Every 30 seconds, the device will remind you to change the cleaning zone by briefly stopping the vibration. Change the cleaning zone after each reminder. The device will turn off automatically after 2 minutes of operation.

Operation modes of the device (see fig. A)

4 – “Soft” mode for delicate, soft brushing.

5 – “Clean” mode for daily teeth brushing.

6 – “Massage” mode for gum massage. Recommended for maintaining healthy gums after brushing your teeth.

7 – “White” mode for teeth whitening. Recommended for people who smoke and drink coffee and tea often.

8 – “Smart” mode for operation of the device according to the scenario set in the AENO app.

Maintenance

WARNING! Do not wash the device and its accessories in a dishwasher. Do not use chemical or aggressive detergents, abrasive pastes, acids, solvents, or metal sponges for cleaning.

WARNING! The device and its accessories should be cleaned after turning off the device and disconnecting the charging base from the power supply.

The brush head bristles should be rinsed under running water (see fig. D). To remove dirt on the handle of the device, rinse it under running water or wipe it with a cloth dampened with water. Afterwards, wipe it with a dry cloth (see fig. D). To remove dirt on the charging base, wipe it with a cloth dampened with water. Afterwards, wipe it with a dry cloth.

Connecting to the AENO mobile app

- Download the AENO mobile app from Google Play or App Store.
- Create a new account or log in to an existing account.
- Turn on the device by pressing the power button (3).
- Make sure that the device battery is not low. The indicator (9) should glow green.

- Put your smartphone as close to the device as possible.
 - On the app home screen, tap “+” in the upper right corner, then tap “Add Device”.
- NOTE.** If the system prompts you to turn on any features on your smartphone, follow the prompts.
- The device will be found automatically, and its image will appear in the “Autosearch” section. Tap the “Add” button.
 - In the next window, tap “+” next to the image of the toothbrush to add it to your account.
 - Wait until the message “Successfully added” appears.
 - Tap “Done”. You can now control the device from the app.

For more information on the device control options in the app, refer to the full manual at aeno.com/documents.

Troubleshooting

- The toothbrush is not turning on. Possible cause: the battery is discharged. Solution: charge the device (charging time – at least 5.5 hours).
- Short operating time after a full charge. Possible cause: the device has not been used for more than 6 months and the battery may not hold a charge well. Solution: contact a service centre to replace the battery.

WARNING! If none of the possible solutions have helped you eliminate the problem, contact the supplier or a service centre. Do not disassemble or attempt to repair the device yourself.

ASBISc reserves the right to modify the device and to make changes and additions to this document without prior notice to users.

The warranty period and service life is 2 years from the date of retail sale of the product. Manufacturer's information: ASBISc Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Made in China. All trademarks and names herein are the property of their respective owners. Up-to-date information and a detailed description of the device, as well as connection instructions, certificates, information about companies that accept quality claims and warranties, are available for download at aeno.com/documents.

SIMPLIFIED DECLARATION OF CONFORMITY. Hereby, ASBISc declares that the device is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address:

aeno.com/userfiles/files/Certificates/CE/DBIS-DB2S.pdf.

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ENG Recycling information These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or turned in to a local recycling facility. You should contact your local household waste disposal service for details.

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