aquapick[®]

CORDLESS ORAL IRRIGATOR

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User Guide

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Advanced Oral Irrigation System
Sistem avansat de irigare orală

AQ-205

AQ-205 is a rechargeable Oral Irrigator.
Please notice the following issues before use.

DANGERS

- ❶ If you do not understand the operating procedures and safety instructions do not use this product to avoid injury caused by improper operation.
- ❷ Please do not operate the product in any way other than this manual.
- ❸ Do not disassemble the product or change the parts of the product to avoid damage to the product or high voltage electric shock. These acts will make the warranty invalid.

WARNINGS

- ❶ This product is only used for teeth cleaning and gum massage. Please do not use it for other purposes.
- ❷ Do not intrude the product into the water to avoid failure, fire, electric shock, explosion and other hazards.
- ❸ If the USB cable falls into the water, dry it before use.
- ❹ Do not plug in or pull off the power plug with wet hands.
- ❺ Do not damage or modify the power cord. Do not stretch, bend or twist the power cord. Do not put weight on the power cord.
- ❻ Please unplug the USB line before maintenance.
- ❼ Do not place the product in a hot or humid place.
- ❽ Please consult the dentist for further use if there is a massive bleeding in the gums after use or if there is still bleeding or other discomfort after one week of use.

ATTENTIONS

- ❶ This product is not suitable for children, infants and people who do not know how to operate unless they are supervised or instructed by their guardian in the use of this product. Please keep it out the reach of children and infants. It cannot be used as a toy.
- ❷ Please contact the manufacturer or the dealers if you need to repair.
- ❸ If you have had oral surgery in the last two months, consult your doctor before using it.
- ❹ This product meets the safety standard of battery device. Involving in pacemaker or other implanted device, consult your doctor before use.

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
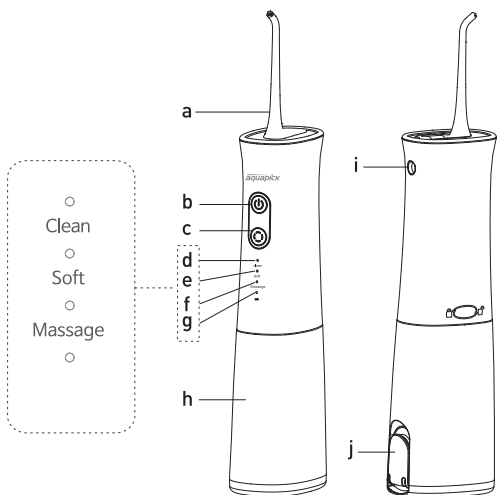
Note:


There may be some liquid remaining inside the main appliance before use (tank, water suction tube, etc). This is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

HOW TO CHARGE

Below devices can be used to charge the oral irrigator (Fig.1)

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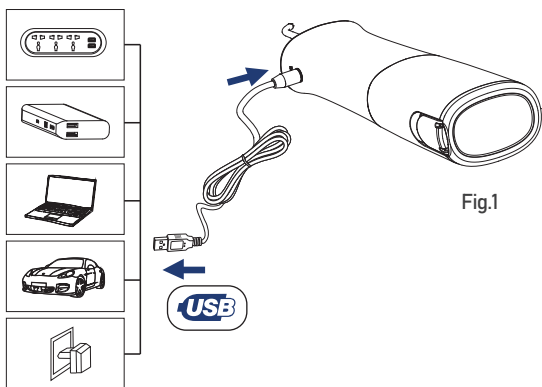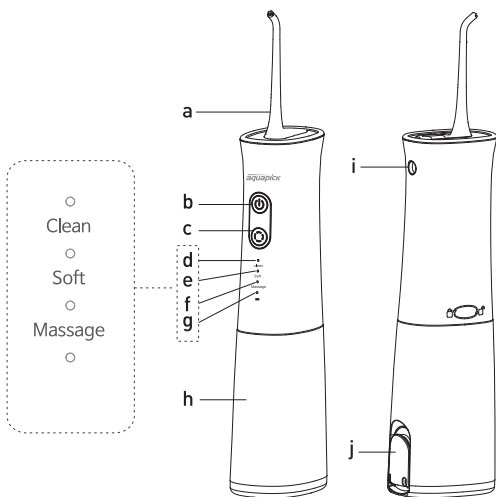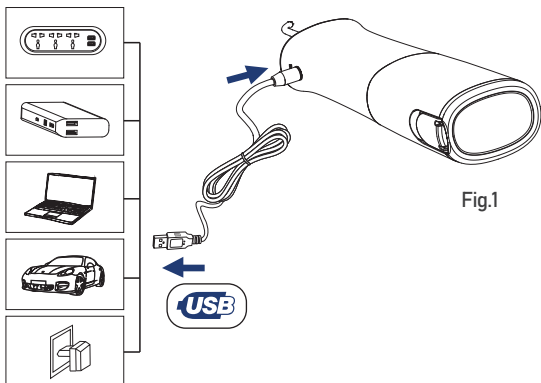
Note:


- When charging the product, please lay the unit on the horizontal table, not standing, to avoid the risk of the unit fall. (as shown in Figure 1)
- While the product is being charged, the indicator is lit in green and keeps blinking. The constant green LED, stops blinking and stays on, indicating that the battery is fully charged.
- If the indicating is blinking green during use, it means the battery is very low. When the battery runs out, the mode light and the motor turn off at the same time, the green light flashes for 10 seconds and then shuts down.

Please recharge the appliance for continuous use.


- It takes 4 hours until fully charged.
- Make sure there is no liquid remaining in the electrical outlet before charging.

USE THE ORAL IRRIGATOR



- 1 Open the water tank cap to fill the tank with water which is not higher than 40°C, then close the cap. (Fig. 2)

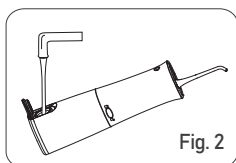- 2 When you use it for the first time, please add small amount of water to water outlet. (Fig. 3)

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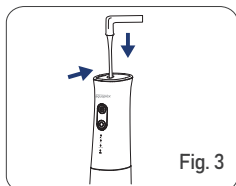
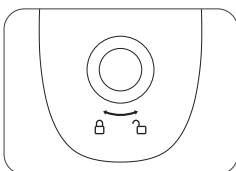
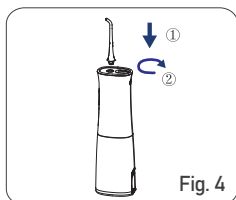
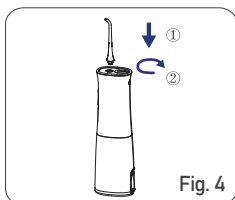
Fig. 4

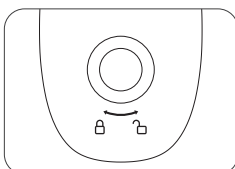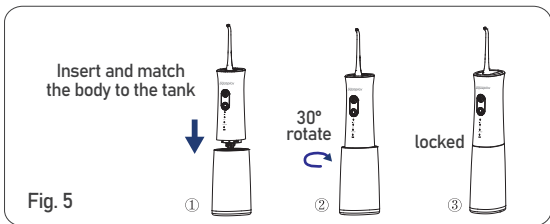④ How to put on the tank. (Fig. 5)

Insert and match the body to the tank; rotate 30° and the tank is locked.

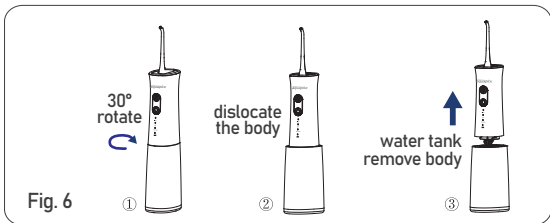
⑤ How to remove the tank. (Fig. 6)

Rotate 30° and dislocate the body and water tank; hold the tank still and pull up the body

⑥ How to operate (Fig. 7)

Press the mode button to switch the mode in order of Clean → Soft → Massage.

After selecting the mode you like, press the On / Off button to start working. Press the On / Off button again to stop working.

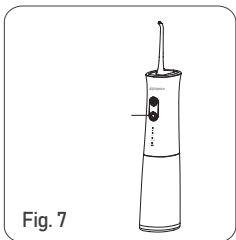Press the On / Off button to start working; Press the mode button to choose the cleaning way you like.


Clean Mode: Deeply clean teeth and massage gums with high pressure.

Soft Mode: Gently clean teeth.

Massage Mode: Efficiently remove stubborn dental plaque and food debris.

⑦ Keep the tip of the jet tip straight toward to your teeth

⑧ Keep your mouth open slightly so that the water can flow out.

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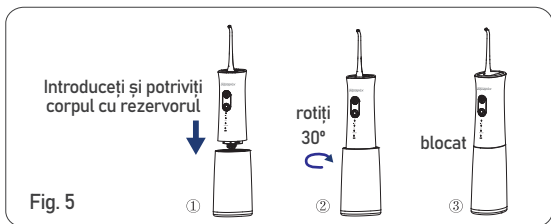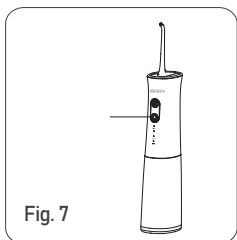
9 Press "On / Off" button to start irrigating (Fig. 8)

- 1) Directing the water jet flow at a 90-degree angle to the teeth and gums.
- 2) Directing the water jet softly along the teeth.
- 3) Directing the water jet along the gum line with the jet tip at a 90-degree angle to the gum tissue.
- 4) Cleaning between the teeth and around any braces, crowns or bridgework thoroughly.

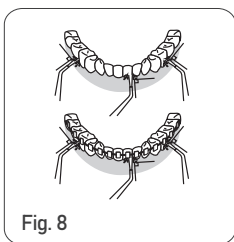
Fig. 8

CLEANING GUIDE

- 1 Press "On / Off" button to turn off the unit.
- 2 Rotate the jet tip 45° counterclockwise and detach the jet tip from the unit. (Fig. 9)
- 3 Open the water tank cap and dispose of the water in the tank.
- 4 Press "On / Off" button to turn on the unit, and drain the water off completely.
- 5 Wipe the device with dry cloth. (Fig. 10)

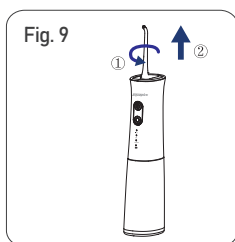
Fig. 9


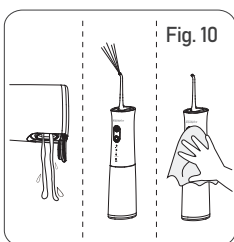
Fig. 10

Note:

To prevent bacteria from propagating, make sure there is no remaining water inside the unit. If you're not going to use it for a long time, be sure to wipe and dry the appliance ASAP before storing it.

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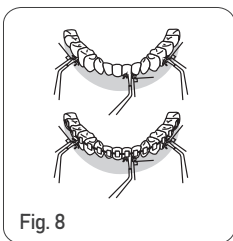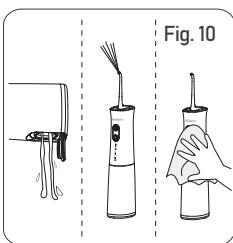
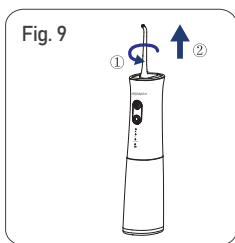
MAINTENANCE

- Use only water or neutral cleaning agent to clean the device.
- Never use caustic or abrasive agents (e.g: vinegar or lime scale removers), since they may damage the device.
- Do not wash by hot water over 40°C.
- Keep the device away from high temperature area and avoid direct sunlight.
- Replace the jet tip every 6 months.

TROUBLE SHOOTING

| PROBLEM | POSSIBLE CAUSE | SOLUTION |
|---|--|---|
| It does not work | You've just purchased the product or it has been left unused for more than 3 months. | Charge the unit for about 4 hours. |
| It only works for a few minutes even after charging | The battery has reached the end of its life cycle. | Contact service center for support. |
| | Charging time was insufficient. | Charge the unit for about 4 hours. |
| Water pressure is very low | There is dirt stuck in the jet tip. | Replace the jet tip. |
| | The jet tip is deformed. | |
| | There is dirt stuck in the filter. | Clean the filter or purchase the filter separately. |
| Water does not come out | The tank is empty. | Fill up tank. |
| | The product is tilted too much while in use. | Use the product in a vertical state. |

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