How to measure body temperature

- Press the on/off button.
- After 5 seconds put the thermometer on the place you want to take temperature - into mouth, under armpit or inside rectum. 5 SECONDS
- When the temperature measuring is complete 10 beeps will sound. IF THE MEASUREMENT IS UNDER THE ARMPIT THE BEEP SIGNAL IS FOR THE APPROXIMATE TEMPERATURE – KEEP THE THERMOMETER UNDER THE ARMPIT FOR 1_2 MINUTES .

Don't take the thermometer out from your body before you hear sound for the end of the measurement, otherwise the LCD will display [- - - °C] and short signals will sound, to warn unfinished or incorrect measurement. Turn the thermometer off, and then turn it on again for new measurement.

Instructions for correct measurement

ZANO digital thermometer is suitable for measurement of body temperature in the mouth, under the armpit and into the rectum. The results on the different spots could vary. The difference* is approximately half a degree.

* The normal temperature and the variation of the body temperatures in the different spots is individual. Measure your body temperature for at least 2 weeks to define yours, on the same spot and time.

The thermometer should be kept at a temperature above 20°C . Precise measurement can be done even if the thermometer is kept on a temperature below 20°C , but in that case the thermometer should be kept on the measurement spot at least 1 min after the signal.


Methods to measure the temperature

Oral:

Measuring time: approximately 10 sec.
(Appropriate for age four or older)

- Place the thermometer tip well under the tongue as indicated by
- Hold the thermometer in the same spot with your mouth closed until the signal sounds.

Do not drink hot or cold fluids, smoke, or exercise 5 minutes before the reading.

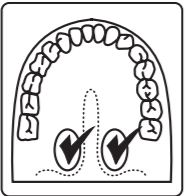
Underarm (Axillary)

Place the thermometer tip in the armpit. The LCD screen needs to face the side as indicated in the drawing.

Firmly fold the arm to hold the thermometer in place and make sure the tip is in contact with the skin.

Before measurement, make sure the armpit is dry.

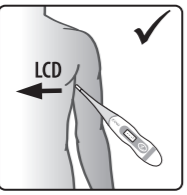Rectal

Measuring time: approximately 10 sec.

- Make sure the thermometer is clean.
- Lubricate the probe with water soluble jelly.
- Gently slide the probe about 1 cm into rectum.

DO NOT FORCE THE TIP INTO THE RECTUM IF RESISTANCE IS ENCOUNTERED.

Indications on the display

- Initial self testing after the switch on.
- Signal for the beginning of the measurement. If the room temperature is above 32°C instead of Lo, the room temperature will be indicated.
- The temperature is more than 42°C.
- Measurement is shown in Celsius degrees.
- Measurement is shown in Fahrenheit degrees.
- (blinking) The thermometer is measuring.
- Memory function. Shows the readout of the last measurement.
- Low battery. Replacement is needed.
- The temperature measurement is incorrect. Turn the thermometer off, and then turn it on again for new measurement.
- The unit does not work properly. Please contact your retailer.

How to switch from Fahrenheit to Celsius


In power off mode, press and hold the on/off button until the display indicates °C or °F

How to replace the battery

When the battery is exhausted indication appears on the display.

To replace the battery:

- Slide back the battery cover.
- Carefully remove the battery from its holder.
- Insert a new battery (LR41)
- When replacing the battery cover, make sure that the seal is not damaged

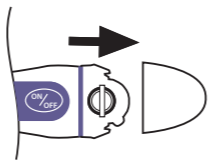Cleaning and disinfection

The thermometer can be cleaned using medical alcohol (95%). Wipe it with clean dry cloth or dry medical cotton before you use the product.

Technical specifications

Measurement range:	32°C – 43°C
Accuracy:	± 0.2°C less than 35.5 °C ± 0.1°C 35.5 °C – 42.0°C ± 0.2°C 42.0°C – 43.0°C
Display unit:	0.1°C
Battery:	One 1.55 V button size battery (LR41)
Power consumption:	0.15 MW in measurement mode
Battery life:	Approx. 1200 measurements per year
Weight:	Approx. 12.6 grams including battery
Dimensions:	129mm x 19mm x 11mm
Storage temperature:	-20°C – 55°C
Operation temperature:	5°C – 40°C
The thermometer is waterproofed.	

Warnings:


The thermometer should be stored at room temperature, in a dry place, away from sources of heat and out of direct sunlight. Strong electromagnetic fields could prevent correct functioning. Keep the thermometer out of reach of children. Never leave children alone while taking a temperature measurement. If temperature rises, consult your doctor. Do not swallow the battery. Do not bite or pull the unit. Do not recharge the battery. Do not dispose the battery in fire. Do not drop the unit. Do not boil the sensor. Do not attempt to disassemble the unit except battery replacement. After eating, drinking or exercising, do not measure the temperature. We recommend verifying the accuracy by authorized laboratory every 2 years.

Applying standards:


93/42/EEC MDD, IEC60601-1:1990, IEC 60601 – 1 – 2:2001, EN 12470 – 3:2000, ASTM E1112:98

Guarantee

As long as the instructions are followed, your Digital Thermometer should last for long years. In the event your Zano Digital Thermometer fails to operate due to manufacturing, materials or workmanship related defects, we will provide you with a replacement. The warranty does not cover damage caused by improper use. Any claim for damages covered by the warranty must be accompanied by the purchase receipt (cash register receipt).

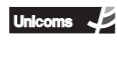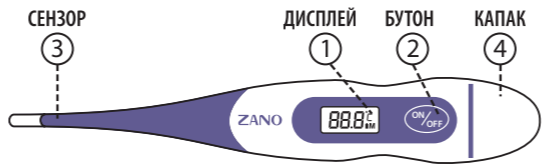Kangfu Medical Equipment Factory, 380 Ningtang East Road, 325600 Yueqing, Zhejiang, People's Republic of China
Luxus Lebenswelt GmbH, Kochstr. 1, Willich, Germany

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